

Agilium® Vantage

Less Pain. More Life.

For 100 years, Ottobock has developed high-quality medical devices to enable patients to regain mobility and independence. Our portfolio of off-the-shelf bracing and supports is designed to improve functional outcomes for patients so they can resume their daily activities and active lifestyles with confidence.

Osteoarthritis (OA) is one of the most common joint disorders worldwide.¹ As the aging population dramatically increases, more people are facing mobility challenges. When it affects the knee, osteoarthritis is one of the most frequent causes of difficulty with walking and climbing stairs.

The utilization of a soft knee brace can be one component in the treatment of osteoarthritis to reduce pain, improve ambulation, and increase stability in patients with knee OA. Drawing from a decade of experience, Ottobock developed the Agilium Vantage, a soft knee osteoarthritis brace for everyday use, recreational activities, and sports. This addition to the Agilium portfolio of knee orthoses dynamically unloads the knee with a single strap so your patients can feel less pain, and experience more life.

¹Grazio 2009

Ottobock US
P 800 328 4058
F 800 962 2549
professionals.ottobockus.com

Ottobock Canada
P 800 665 3327
F 800 463 3659
professionals.ottobock.ca

©2019 Otto Bock HealthCare LP · 16454 · 09/19

ottobock.

Agilium® Vantage

Less Pain. More Life.



Quality for life


Information for practitioners

Overview

Meet the Agilium Vantage, a low-profile wraparound knee brace designed to provide lasting pain relief for patients with mild to moderate unicompartmental osteoarthritis. The dynamic Y force strap system unloads the affected compartment of the knee joint with one pull, while reducing rotation and providing stability.

Indications

- Mild to moderate unicompartmental knee osteoarthritis
- Medial or lateral knee pain
- Rheumatoid knee arthritis
- Post-operative meniscus reconstruction or ligament injuries needing unilateral relief

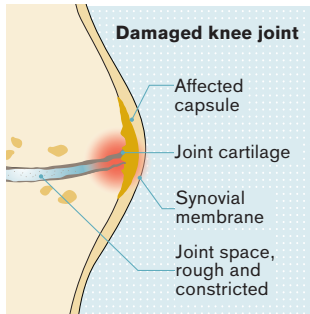

Knee osteoarthritis pain

- Bone spurs
- Cartilage loss
- Bone exposure
- Meniscus erosion
- Narrowing joint space

Treatment with the Agilium Vantage

- Utilizes three-point principle to stabilize the knee and provide pain relief
- Single upright design with four anterior calf and thigh wraps

Dynamic strap system advantage

- Dynamic Y force strap unloads the affected medial or lateral compartment of the knee
- Reduces rotational forces

Features & Benefits



Dynamic Y force strap system

Unloads the knee joint with a single pull to minimize rotation and reduce pain




Low-profile design

Single metal upright with hinge for fixed flexion and extension limits



Easy to don and doff

Four adjustable, numbered closure straps for quick and easy application

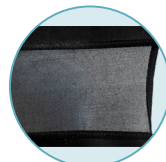

Breathable material

Lightweight material prevents popliteal bunching for user comfort and increased compliance



Reduce inventory

Trimmable sleeves to fit a broader range of patients with three sizes

Ordering



How to Order

1. Order online

US: professionals.ottobockus.com

CA: professionals.ottobock.ca

2. Order through Ottobock Customer Service

US: 800 328 4058 or USCustomerService@ottobock.com

CA: 800 665 3327 or CACustomerService@ottobock.com

3. Order through your Ottobock Sales Rep

Standard Warranty

6 months

Get in Touch

Don't hesitate to contact us if you have further questions or would like more information.

Sizing						
	Article Number	Side	Size	A (in)	B (in)	Type
A	50K306	L/R	XS-S	14-18	11-14	MO/LO
B		L/R	M-L	18-23.5	14-18	MO/LO
		L/R	XL-3XL	23.5-32	18-24	MO/LO

Example format: 50K306=L-XS-S-MO or 50K306=R-M-L-LO

Place hinge on the side of the affected compartment