

Agilium® Softfit

Less Pain. More Life.

Ottobock's Agilium portfolio is designed to support patients with knee osteoarthritis (OA) for everyday use, recreational activities, and sports. The Agilium Softfit (50K90) is a single upright knee orthosis with a textile base body for patients with unicompartmental knee OA.

Features & Benefits

- Proven three-point pressure system causes direct knee offloading and reduces pain
- Adjustable bars for quick and straightforward fitting to treat both varus and valgus types
- Comfortable, slim design and patented vector grip material ensures firm fit without slipping
- Easy to don and doff by the patient themselves


Indications

Medial or lateral unicompartmental knee osteoarthritis.



Adjustments

The corrective force can be set with an Allen wrench above the joint to easily adjust the varus/valgus angle of the orthosis.

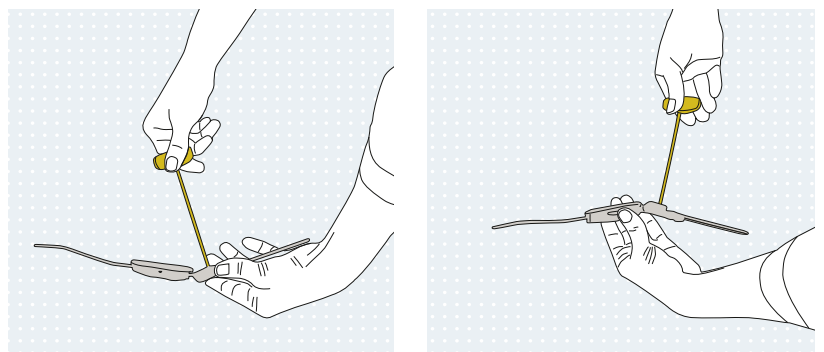

Ordering Information

Article number	Side	Size	Circumference	
			A (cm)	B (cm)
50K90	L/R	S	35-38	44-48
	L/R	M	38-41	48-52
	L/R	L	41-44	52-56
	L/R	XL	44-48	56-61
	L/R	XXL	48-51	61-67

Example format: 50K90=L-S or 50K90=R-M