

ottobock.

Agilium[®] Forte

Less Pain. More Life.






Quality for life

Overview

Meet the Agilium Forte, a single upright rigid knee orthosis designed to provide lasting pain relief for patients with moderate to severe unicompartmental osteoarthritis.

Indications

- Moderate to severe unicompartmental knee osteoarthritis
- Medial or lateral knee pain
- Rheumatoid knee arthritis
- Post-operative meniscus reconstruction or ligament injuries needing unilateral relief

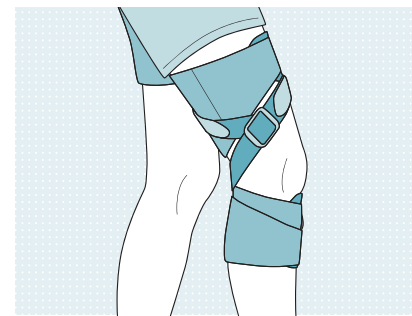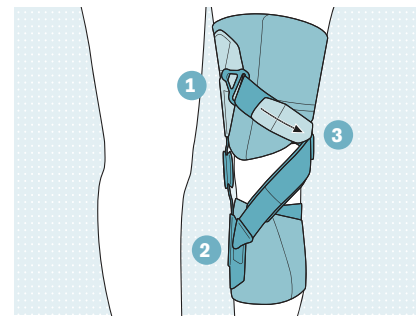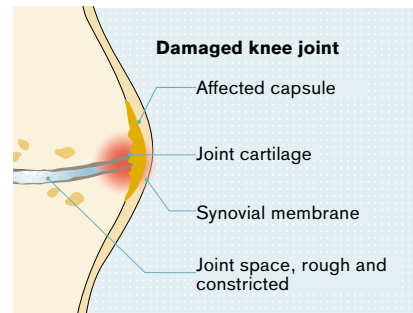

Knee osteoarthritis pain

- Bone spurs
- Cartilage loss
- Bone exposure
- Meniscus erosion
- Narrowing joint space

Treatment with the Agilium Forte

- Utilizes the three-point principle to stabilize the knee and provide pain relief
- Single upright design with rigid shells and three anterior calf and thigh wraps

Dynamic strap system advantage

- Patented Y force strap unloads the affected medial or lateral compartment of the knee
- Anti-rotation strap reduces orthosis migration

Less Pain. More Life.


For 100 years, Ottobock has developed high-quality medical devices to enable patients to regain mobility and independence. Our portfolio of off-the-shelf bracing and supports is designed to improve functional outcomes for patients so they can resume their daily activities and active lifestyles with confidence.

Osteoarthritis (OA) is one of the most common joint disorders worldwide.¹ As the aging population dramatically increases, more people are facing mobility challenges. When it affects the knee, osteoarthritis is one of the most frequent causes of difficulty with walking and climbing stairs.

The utilization of a knee orthosis can be one component in the treatment of osteoarthritis to reduce pain, improve ambulation, and increase stability in patients with knee OA. Ottobock developed the Agilium Forte, a rigid knee osteoarthritis brace for everyday use, recreational activities, and sports. This latest addition to the Agilium portfolio of knee orthoses dynamically unloads the knee with a single pull so your patience can feel less pain, and experience more life.

¹Grazio 2009

Features & Benefits

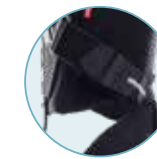


Dynamic Y force strap system
Patented dynamic Y force strap unloads the knee with one pull to provide pain relief while the anti-rotation strap reduces orthosis migration on the leg



Low-profile design
Lightweight, rigid shells for stability and open sleeve design to reduce popliteal bunching for user comfort and increased compliance



Highly adaptable
Single upright with three adjustable height settings to fit a broad range of patients

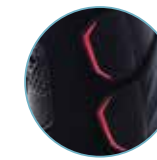

Easy to don and doff
Three adjustable, numbered closure straps for quick and easy application



Reduce inventory
Trimmable sleeves reduce inventory to three sizes

Ordering



How to Order

1. Order online

US: professionals.ottobockus.com
CA: professionals.ottobock.ca

2. Order through Ottobock Customer Service

US: 800 328 4058 or USCustomerService@ottobock.com
CA: 800 665 3327 or CACustomerService@ottobock.com

3. Order through your Ottobock Sales Rep

Standard Warranty

6 months

Get in Touch

Don't hesitate to contact us if you have further questions or would like more information.

		Sizing					
		Article Number	Side	Size	A (in)	B (in)	Type
A	6 in	50K305	L/R	XS-S	14-18	11-14	MO/LO
			L/R	M-L	18-23.5	14-18	MO/LO
B	6 in		L/R	XL-3XL	23.5-32	18-24	MO/LO

